


Thermal analysis of the sample#P9830-EOPOEO

Thermal analysis of the samples was carried out on a TA Q100 differential scanning calorimeter at a heating rate of 10°C/min. The midpoint of the slope change of the heat flow plot of the second heating scan was considered as the glass transition temperature (T_g).

Thermal analysis results at a glance (EO-PO-EO)

Sample	T_m (°C)	T_c (°C)	T_g (°C)
EO block	53	29	-57
PO block		-	-


Typical thermogram for the PO-EO block

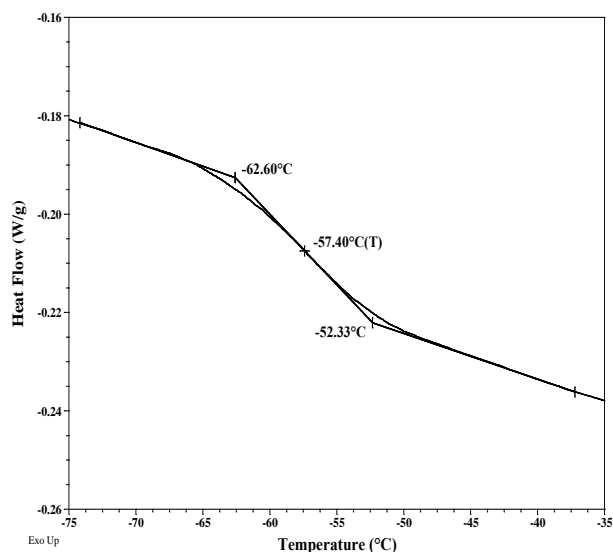

Melting and crystallization curve for the sample:

The melting temperature (T_m) was taken as the maximum of the endothermic peak where as the crystallization temperature (T_c) was considered as the minimum of the exothermic peak.

Melting curve for PEO block:

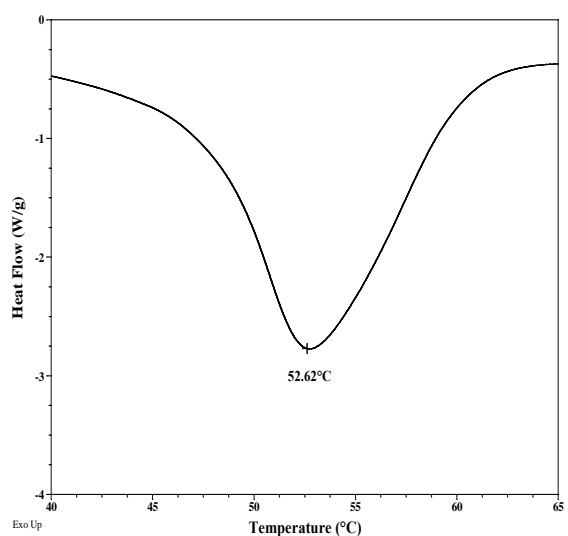

Crystallization curve for PEO block:

