


Sample Name: Polystyrene

Sample #: P8097-S

Structure:

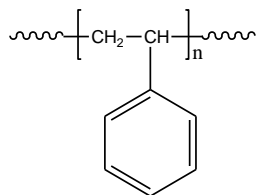


Composition:

Mn x 10 ³	PDI
29.0	1.10

Synthesis Procedure:

Polystyrene is obtained by living anionic polymerization of styrene as illustrated below:

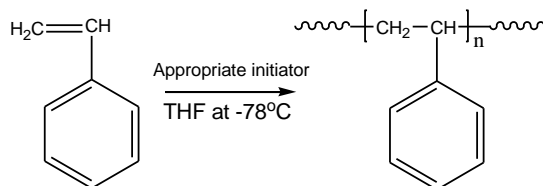

Characterization:

The molecular weight and polydispersity index (PDI) are obtained by size exclusion chromatography (SEC) in THF. SEC analysis was performed on a Varian liquid chromatograph equipped with refractive and UV light scattering detectors. Three SEC columns from Supelco (G6000-4000-2000 HXL) were used with triple detectors from Viscotek Co.

Thermal analysis:


Thermal analysis of the sample was carried out using a differential scanning calorimeter (TA Q100) at a heating rate of 10°C/min. The inflection glass transition temperature (T_g) has been considered.

Solubility:

Polystyrene is soluble in DMF, THF, toluene and CHCl₃. It precipitates from methanol, ethanol, water and hexanes.

SEC of Homopolymer:


P8097-S

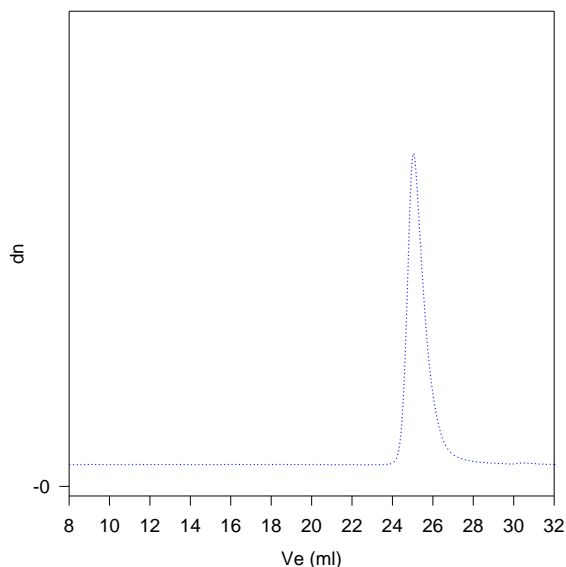

Size Exclusion Chromatography of polystyrene

M_n=29000, M_w=32000, M_w/M_n=1.10

¹H NMR of the Product:

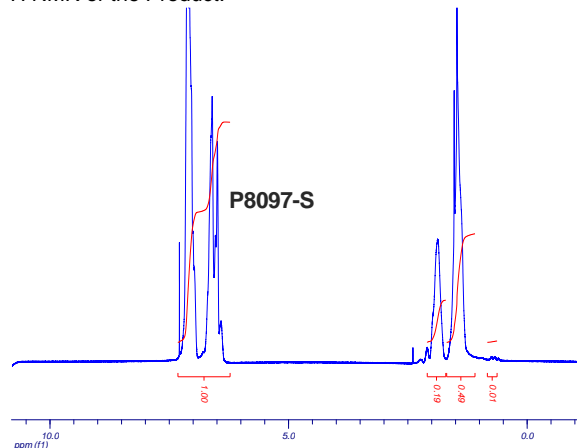

T_g of polystyrene as function of molecular weight

