

Thermal analysis of the sample# P2591-DMStBuA

Thermal analysis of the samples was carried out on a TA Q100 differential scanning calorimeter at a heating rate of 10°C/min. The midpoint of the slope change of the heat flow plot of the second heating scan was considered as the glass transition temperature (T_g).


Melting and crystallization curve for the sample

The melting temperature (T_m) was taken as the maximum of the endothermic peak where as the crystallization temperature (T_c) was considered as the minimum of the exothermic peak.

Thermal analysis results at a glance

Sample	T_m (°C)	T_c (°C)	T_g (°C)
DMS block	-49	-	-127 (Lit.)
tBuA block	-	-	23

Thermogram for tBuA block

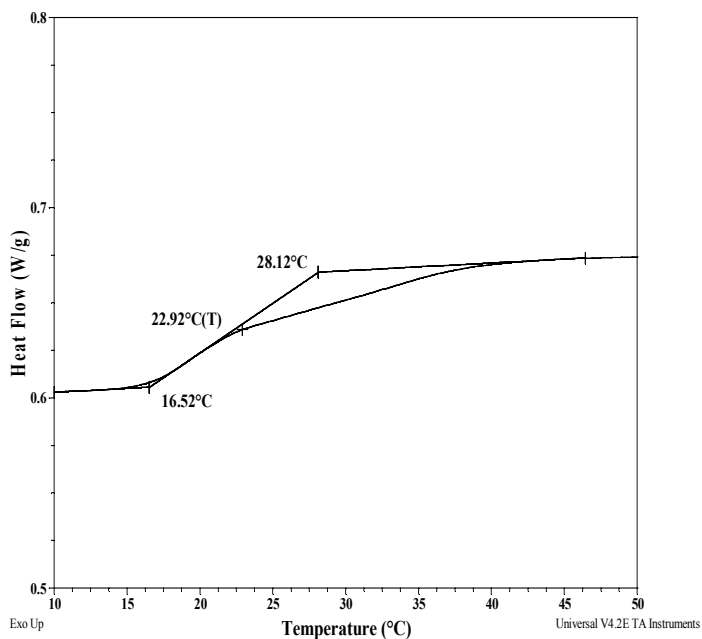

Melting curve for DMS block:

